


► Thank you for choosing **SurgiTel®**!

Your SurgiTel® LED headlight has been designed to meet the demanding needs of today’s clinical professionals by allowing you to see, work, and feel your best.

► TECHNICAL DATA

LED Light	Several optics available
Battery Charge Life (with 3 batteries)	Up to 4 hours
Battery	Rechargeable Lithium Ion Battery, 3.7V 350 mAh
Battery Charger (cradle) & AC Power Adapter	Input: 100-240V, 50-60Hz, 0.35A Output: 5V  200mA

► GENERAL WARNINGS

- The user of this product should be thoroughly familiar with the use and care of this product.
- The user should carefully study the manual before making any attempt to use this product clinically.
- Before each procedure carefully inspect the LED system to ensure proper maintenance.
- Use care not to point the light directly at the eyes during operation; patient use of safety glasses is recommended.

FIRE HAZARD: DO NOT PLACE THE LED LIGHT UNIT ON ANY FLAMMABLE OBJECT WHILE IT IS OPERATING.

WARNING:

- The LED light unit will get hot when left on. Handle with care. It is recommended to turn off the LED before handling LED light or making any adjustments to it.
- Avoid direct light contact with the eyes.
- Avoid direct light irradiation in patient’s eyes; patient use of safety glasses is recommended.

Battery Charger (cradle) & AC Power Adapter

- Use only the Battery Charger & AC Power Adapter supplied with this system
- Do not disassemble, modify or short circuit
- Do not leave in hot places, such as inside a car or outside in hot weather
- Do not allow to get wet



► WARRANTY

SurgiTel’s mission is to provide the best products for your *Vision & Ergonomics at Work*. To support this commitment SurgiTel® warranties all of its products.

SurgiTel® strives to make the best loupes and light systems available. Please inspect and try on product immediately upon receipt. If you notice any issues, contact our customer service department at 800-959-0153.

For both warranty and out-of-warranty issues, contact SurgiTel® Customer Service regarding repairs.

Warranties vary by product. Check www.surgitel.com/warranty for more details. Check www.surgitel.com/help for up to date care information/instructions. Please keep up to date on care information/instructions, as warranties can be affected if product is misused.

► PATENT


True color beam family optics, U.S. Patent Numbers: US 9,383,077 B2; US 9,494,299 B2

SurgiTel®
a division of General Scientific Corporation
77 Enterprise Drive • Ann Arbor, MI 48103-9503
E-mail info@surgitel.com • www.surgitel.com







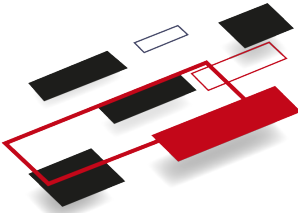

Odyssey Wireless LED Air

USER INSTRUCTIONS



➤ CLASSIFICATION

- Protection against electric shock: Class II including AC/DC Adapter.
- Applied parts : No applied parts.
- Protection against harmful ingress of water: ordinary IPXO.

	Direct current		Class II equipment
	CE mark		Serial No.
	Alternating current		Attention, consult accompanying documents

➤ MAINTENANCE

Storing the LED system in a safe place will prolong the life of the unit. The following guidelines will help in sustaining the LED systems useful life:


- Keep the optical face on the light unit from touching hard surfaces which may cause scratches on the surface. Scratches on the surface will diminish the light quality.
- When cuts or punctures are visible on the system’s light cord, the wireless LED system is no longer safe to use; it should be taken out of service immediately.

Fixing Hinge Tension (too loose or too tight):

Step 1: Use provided hex key to remove hinge set screw (see photo to the right)

Step 2: Use two flat-head screwdrivers (placed on each side) to loosen/tighten special screw and nut

Step 3: Once desired tension is met, reassemble and tighten set screw

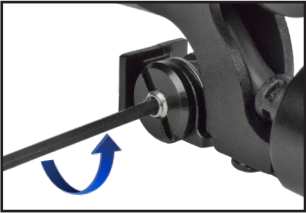

➤ CLEANING




IMPORTANT

- Do not immerse in water or any liquid; wipe down with damp cloth. Do not use HOT or COLD water; using “room-temperature” water is acceptable
- Avoid using running water and disinfectant sprays on or near open slots located on light unit. Carefully wiping with damp cloth is ok.
- Do not use Lysol® or any disinfectants containing glutaraldehydes, iodophors, ammonia, alcohol, or ethanol
- Do not autoclave

Disinfection: Use SciCan® “OPTIM 33 TB” wipes/solution only.
*Note: SciCan® OPTIM 33 TB solution can be purchased through SurgiTel in 8 oz. bottles

FACTORY REFURBISHMENT CLEANING: is available if you come across a difficult cleaning situation. Call customer support for details: 800-959-0153

➤ OPERATION

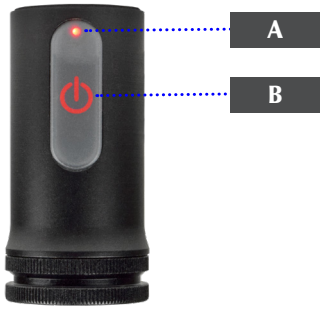

HOW TO TURN LIGHT ON AND OFF

Touch on/off symbol (B) on control panel to turn light on or off. When light unit is turned on, a red indicator light (A) will appear.

BRIGHTNESS OF LIGHT

Brightness of light is nonadjustable and is factory set at max brightness.

Note about light-cured composites

- Using the appropriate COLOR FILTER (orange filter) is STRONGLY RECOMMENDED when using light-cured materials with the LED Light

➤ OPERATING CONDITION

- Ambient temperature range 10 to 40°C
- Relative humidity range 30 to 75%
- Atmospheric pressure range 70 to 106 kPa

➤ TRANSPORT & STORAGE CONDITION

- Ambient temperature range -40 to 70°C
- Relative humidity range 10 to 100%
- Atmospheric pressure range 50 to 106 kPa

➤ ATTACHING LIGHT ADAPTER TO LOUPE




Dove Tail Adapter



Clip Adapter (provided w/ kit)

- Attach the LED light to loupes or glasses using appropriate adapter

➤ ADJUSTMENTS

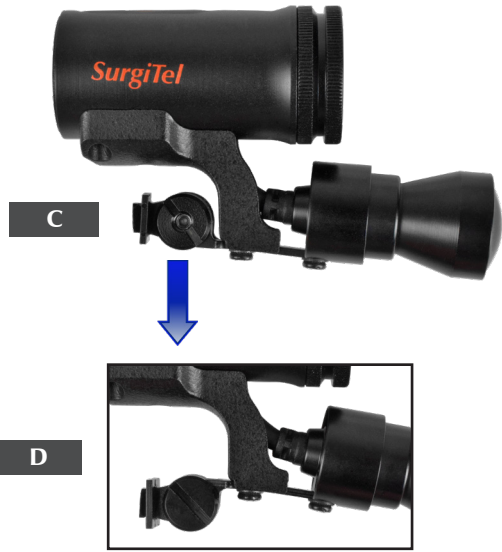

Vertical Alignment (C):

Maximum range of adjustment for light system placement

Optional Light Adapter Assembly (D):

Recommended option if experiencing attachment difficulties or alignment issues with loupe, remove screws and reassemble by changing the light adapter position (as shown)

➤ RECHARGEABLE BATTERIES



Charge Level Indicator

- Charge Time = approximately 1 hour
- Charge Level Indicator Lights (3 qty., color: blue): All 3 lights on = fully charged
- Maximum continuous operation (1 battery) from full power is approximately 1.20 hours

Note: To maximize use of Wireless LED Air, remember to replace and charge battery after use